


One Step Strip Style Alcohol Saliva Tests

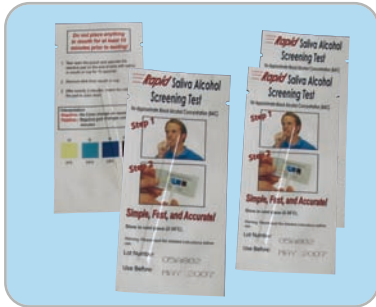

Rapid Alcohol Test is intended as a rapid method to detect the presence of alcohol in saliva for blood alcohol concentration (BAC) greater than 0.02%. Rapid Alcohol test is intended for the semi-quantization of ethyl alcohol in human saliva by comparing the color developed at the end of the reaction period with the color chart provided.

- *Fast and easy to use, no other reagent is needed;*
- *Results are read visually without the need for instrumentation;*
- *Highly sensitive.*
- *The result can be observed within 20 seconds.*

Catalog No.	Product Name	Package	Performance Characteristics	Reading time	Storage Condition	Shelf Time
ALC 711	Saliva Alcohol Test Strip	1 Strip / Pouch 100 pouches / Box	sensitivity: 0.02% alcohol concentration in saliva;	Read result within 2 minutes.	Room temperature (4 to 30°C)	24 Months

One Step Troponin I Rapid test



One-Step Troponin I test is a fast in vitro immunochromatographic test. It is designed for qualitative determination of cardiac troponin I (cTnI) in human serum or plasma specimens as an aid in the diagnosis of myocardial infarction. Principle: sample mobilizes gold anti-cTnI conjugate that is coated on the conjugate pad. The mixture moves along the membrane by capillary action and reacts with anti-cTnI antibody that is coated on the test region.

- *Highly sensitive and specific;*
- *Assays are easy to perform;*
- *Qualitative results are read visually without any instrument;*
- *No other reagent is needed;*
- *Takes less than 15 minutes to be sure.*

Catalog No.	Product Name	Package	Performance Characteristics	Reading time	Storage Condition	Shelf Time
cTnI 721	Troponin-I Serum/plasma Cassette	1 Test/Pack 25 Packs/Outer	Sensitivity: 1.0ng/ml	Read at 15 minutes. Do not read after 15 minutes.	Room temperature (4 to 30°C)	24 Months



Atlas Link Biotech. Co., Ltd

Address: Room 811 Zeyang Plaza, No.166 Fushi Road, Shijingshan Dist. Beijing 100043, China

Tel: 86-10-8890.9112/9113 Fax: 86-10-88909115

Website: www.ivdpretest.com

E-mail: rapidtest@ivdpretest.com



ATLAS LINK, INC.

Address: 8807 Sudley Rd #101, Manassas, VA 20110, USA