


One Step Dengue IgG / IgM Antibodies Test

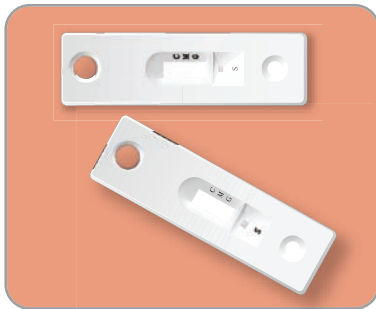

Dengue IgG / IgM Rapid Test is a solid phase immunochromatographic assay for the rapid, qualitative and differential detection of IgG and IgM antibodies to dengue virus in human serum, plasma or whole blood. This test is intended for professional use as aid in the presumptive diagnosis between primary and secondary dengue infection.

- *Highly sensitive and specific;*
- *No other reagent is needed;*
- *Qualitative results are read visually without any instrument.*

Catalog No.	Product Name	Package	Performance Characteristics	Reading time	Storage Condition	Shelf Time
DNG 431	Dengue IgG/IgM Ab Test Cassette(Serum, plasma or whole blood)	25Tests/Box	Sensitivity: 98% Specificity: 100%	Read results within 15 ~ 30 minutes	At room temperature (2-30° C)	24 Months

One Step H. pylori (Helicobacter pylori) Tests



The H. pylori test is a rapid test for the qualitative detection of antibodies of all isotypes (IgG, IgM, IgA, etc) specific to Helicobacter pylori in human serum, plasma or whole blood. This test kit is intended as an aid in the diagnosis of H. pylori infection in patients with gastrointestinal symptoms.

The H. pylori antigen-colloid gold conjugate and serum sample moves along the membrane chromatographically to the test region (T) where H.pylori capture antigen is pre-coated, and forms a visible line as the “antigen-antibody-antigen gold particle complex” forms.

- *Easy to conduct;*
- *qualitative results are read without instrumentation;*
- *Highly sensitive and specific;*
- *Require no additional reagents.*

Catalog No.	Product Name	Package	Performance Characteristics	Reading time	Storage Condition	Shelf Time
H.P 441	Helicobacter Pylori serum/plasma Test Cassette	1 Test/Pack 25 Packs/Outer	Sensitivity: 95.9% Specificity: 89.6% total agreement: 93.4%	Interpret results within 10 minutes; Do not interpret after 20 minutes	Room temperature (4 to 30°C)	24 Months



Atlas Link Biotech. Co., Ltd

Address: Room 811 Zeyang Plaza, No.166 Fushi Road, Shijingshan Dist. Beijing 100043, China

Tel: 86-10-8890.9112/9113 Fax: 86-10-88909115

Website: www.ivdpretest.com

E-mail: rapidtest@ivdpretest.com



ATLAS LINK, INC.

Address: 8807 Sudley Rd #101, Manassas, VA 20110, USA